


Biodex Medical Center

55 red st
Shirley, NY 11967

Daily Calibration w/ Constancy Detail Report

Printed on: 05/08/2017 12:53:46 PM

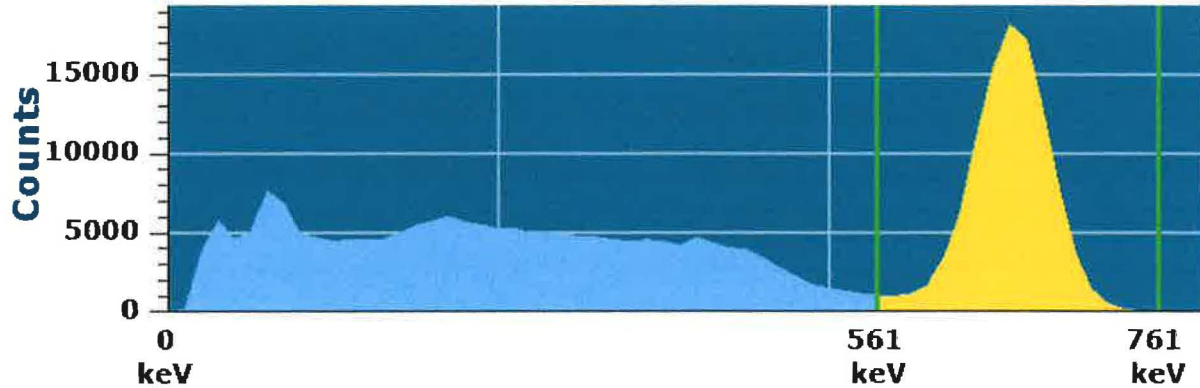

Test Information

Calibration Date/Time : 05/08/2017 9:30:22 AM

Detector S/N : 2007002

Isotope : Cs-137

Staff : C S

Count Time (sec) : 20

662 keV Channel : 51

FWHM : 7.8 %

Constancy Variance : -5.20 %

Acceptable Variance : 5.0 %

Constancy Source Info : Source ID: 1, Activity: 1 μ Ci, Date: 11/29/01

Spectrum Statistics

Total Counts : 284,894

Total CPM : 854,682

ROI Statistics

ROI (keV) : 561 - 761

ROI Counts : 102,508

ROI CPM : 307,524.0

Peak Statistics

Peak (keV) : 662

Peak Counts : 18,114

Peak CPM : 54,342.0

Comments :

Signature : _____

Date Signed : _____