


Wipe Test Detail Report

Printed on: 6/1/2017 3:38:15 PM

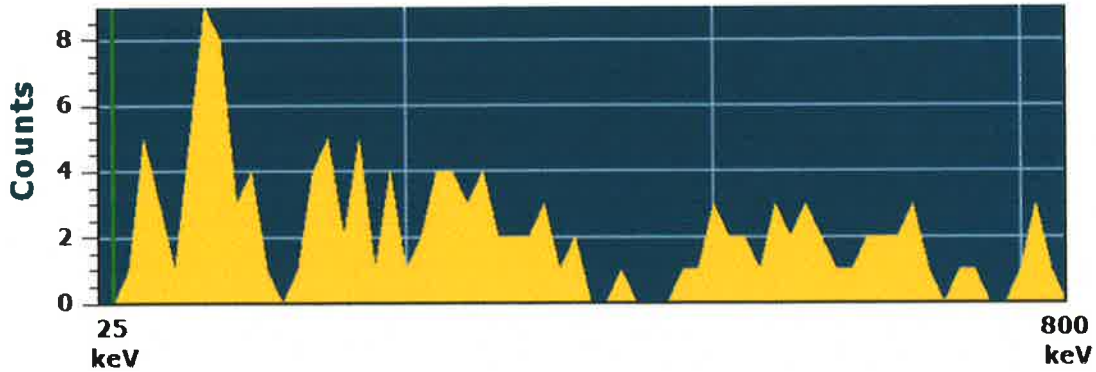

Test Information

Date/Time : 6/1/2017 3:36:32 PM
 Detector S/N : 123
 Staff : Ann Marie Alessi

Full Spectrum

Net Counts : 200
 Total Net CPM : -96.0
 Test Time (secs) : 30
 Net Counts >750 keV : 72

Area Name : **Camera room**
 Location Name : **Camera**
 Status : **PASS - Acceptable Level**

Test Details

Isotope	ROI (keV)	Background CPM	MDA	Net CPM	Net Activity	Trigger	Status
Wide Window	25 - 800	302.00	333 DPM	-44.0	-122 DPM	20,000 DPM	Pass
Co-57	50 - 158	98.00	-7,136 DPM	-30.0	3,034 DPM	20,000 DPM	Pass
F-18	420 - 600	74.00	-6,273 DPM	-32.0	3,236 DPM	20,000 DPM	Pass
Tc-99m	110 - 163	36.00	-4,541 DPM	-18.0	1,820 DPM	20,000 DPM	Pass

Comments :

Signature : _____

Date Signed : _____